



# RC<sup>®</sup> | POWER BASIC

Distributed in the USA by: Dermavista Medical Aesthetics  
Dermavista.com • 1.800.767.6765

THERMO RADIOFREQUENCY + REJUVENATION



## Natural perfection for the skin.

*A system designed to lift and tighten without surgery, restoring the vitality of the skin.* The vanguard application of both cold as a therapeutic agent in combination with radio frequency yields excellent results. By cooling the tissue maximum RF potency can be achieved.



6MHz Facial and Corporal Heads



INTERCHANGEABLE HAND PIECE

REJUVENATION • CUTANEOUS FLACCIDITY • LIFTING WITHOUT SURGERY  
CELLULITE • CREASE AND LINES OF EXPRESSION

### SPECIFICATIONS

#### EASY TO TRANSPORT

A new RC Power model thought to adapt to all of your spaces

Weight: 10 kg

Dimensions: height 41 cm - width 42 cm - depth 46cm

Display: LCD 128x64

Keyboard: Trackpad

### NEW SYSTEM *ELECTRONIC COOLING*

#### RADIOFREQUENCY

RF Frequency: 6 MHz

RF Output: an interchangeable hand piece

Cooler: Electronic Cooler (8°C)





**RC®**  
POWER BASIC

## EQUIPMENT

### Unprecedented sophistication

A new design that give you the assurance of function, comfort and system performance. Built with top quality materials and respecting the **RC Power®** elegant and sober lines.

The evolution of the **RC Power®** permits that the new line of equipments be easily transportable occupying a reduced space inside a consulting room.



### Electronic Cooling

Thanks to the effect of the cold/hot contrast produced by the utilization of the Electronic Cooling, it is possible to diminish the tolerable sensibility, permitting the use of a more potent Radiofrequency. Therefore getting better results in less time.

### High intensity work

### High potency

### Professional intensive use

### Without consumables, without maintenance

### Radiofrequency at 6MHz

The Radiofrequency therapy, in 6MHz is an excellent tool to achieve skin rejuvenation. Thanks to the inclusion and use of cold as a therapeutic agent in combination of the Radiofrequency, excellent results are permitted. This is because by previously cooling the tissue, we can apply the equipment's maximum potency. The resistance of the tissues, as the Radiofrequency passes, is very low so the sensation that the patient perceives during the RC Power application is a marked thermic effect. The bioelectric effect generated by the Radiofrequency is associated with this detachment of the heat produced in the depth of the tissue.

### Visible Results

The effects that are clearly in evidence, immediately after every application, in the dissimulatin of wrinkles, expression creases and improvement of the patient's skin quality. Thanks to the interchangeable heads, Radiofrequency can be applied at different frequencies, which permits the realization of applications at different levels on the same area, therefore potentiating the final result. The immediate effectds achieved are: Hyperaemia, Fibroblastic stimulation and an increase in temperature in dermic levels.

### Radiofrequency 3D



The heat produced in the tissue by the passage of the RF is an excellent alternative to obtain non invasive reaffirmation of flabby and aged skin.




The controlled heating generated at this frequency reaches a superficial cutaneous temperature between 38-42 ° C, capable of stimulating the existing dermic collagen and producing a transition process in the denominated collagen contraction.



These reactions are used to give more firmness and quality to the skin. The result is smoother skin and an evident reduction of wrinkles and flaccidity, including the effect of the orange skin, characteristic in cellulite.



The cooler protects the epidermis from overheating. It offers protection to the superficial layers of the skin from an excessive increase in temperature and also favors the obtention of the denominated contrast effect.

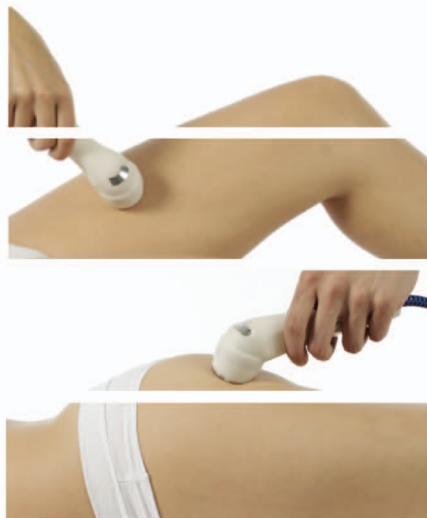

## THERAPIES

